

 **CTLIFT SYSTEM™**  
REINVENTING IMA ACCESS




**TSI** TeDan Surgical  
Innovations

# CT LIFT SYSTEM

## IMA Exposure

TSI is reinventing IMA access for standard and minimally invasive procedures by improving on proven technology with a modern design that provides optimal exposure and more positioning versatility. Cut the cable – the sleek design eliminates use of braided steel cables which simplifies cleaning and reduces required maintenance.

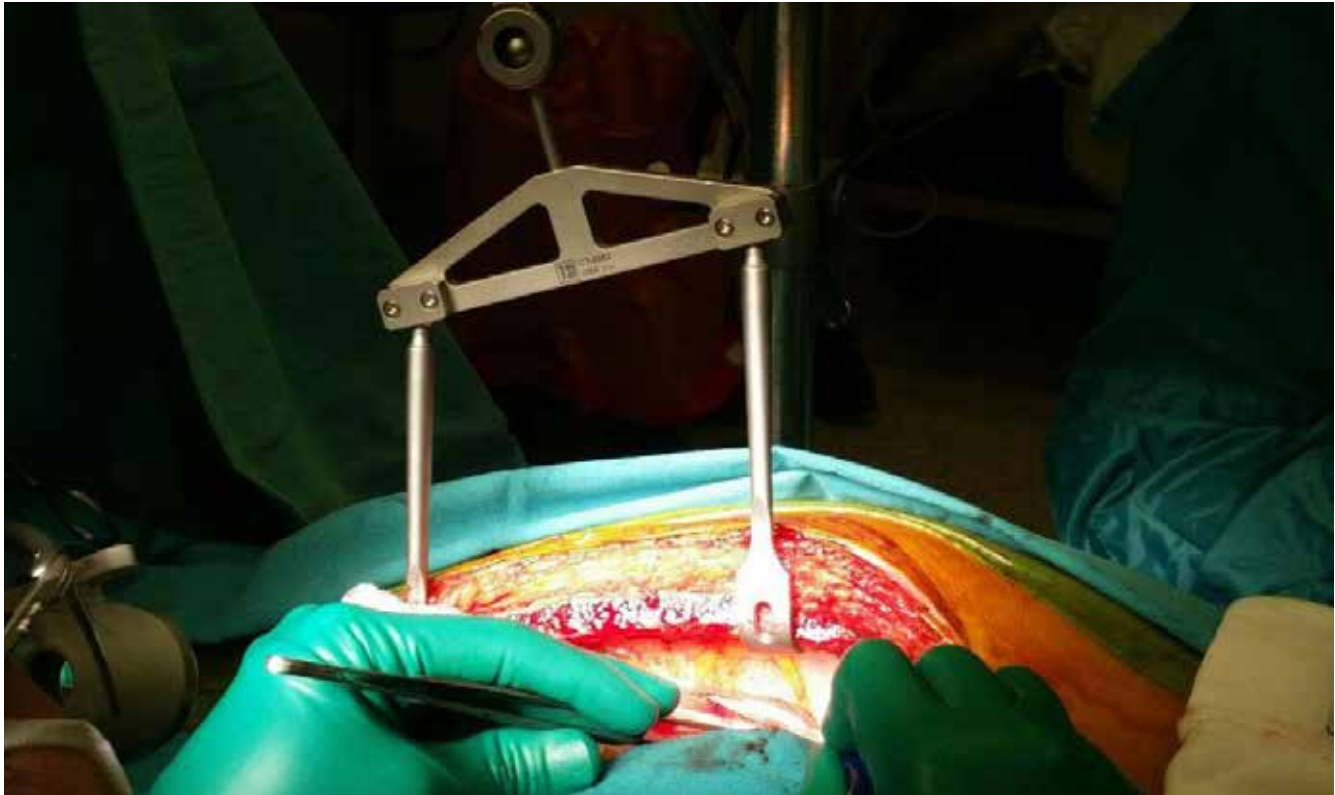

### Sternal Setup Using Dual Hook

#### Open IMA Exposure

CT-0052

#### Dual Hook

- Time tested design to engage the sternum
- Improved connection to lifting mechanism

### MI Thoracotomy Setup Using IMA Blade

#### Minimally Invasive IMA Approach

CT-0050, CT-0300 -CT-0302

#### Blade Hook and IMA Blades

- Wide and narrow IMA blades attach to existing CT-0100 Ultravision™ retractor accommodating diverse patient anatomies
- Pivot function of retractor improves access and visibility

## Versatility Made Simple

Adaptive setup options ensure patient-specific compatibility

Unique adjustment mechanism of extension arm enables various angles and post lengths, providing exact positioning.



## Improving System Assembly

Creating effortless and stable IMA exposure

CT Lift System components simplify the setup process through a unique table locking post and one-step cam mechanisms.



Cam mechanism provides precise positioning along the posts and creates a rock-solid lock in one turn of the handle

Universal hook holder optimizes attachment of preferred IMA approach to the lift mechanism



- Table post attaches directly to OR side rail over the drape
- Tightening handle stows when not in use
- No stainless steel braided cables


## CT-1705 CT LIFT SYSTEM

### Reinventing IMA Access:

- For standard and minimally invasive approaches
- Rapid setup
- Easy to clean
- Lifetime warranty

### Set includes:

Cat. No.	Description
CT-0050	IMA Blade Hook
CT-0052	Open IMA Dual Hook
CT-0300	IMA Blade, Wide
CT-0302	IMA Blade, Narrow
CT-0705	Sterilization Tray, CT-1705
HS-0016	Table Clamp, 450 mm, Purple
HS-0019	Extension Arm
HS-0100	Lift Assembly, Purple

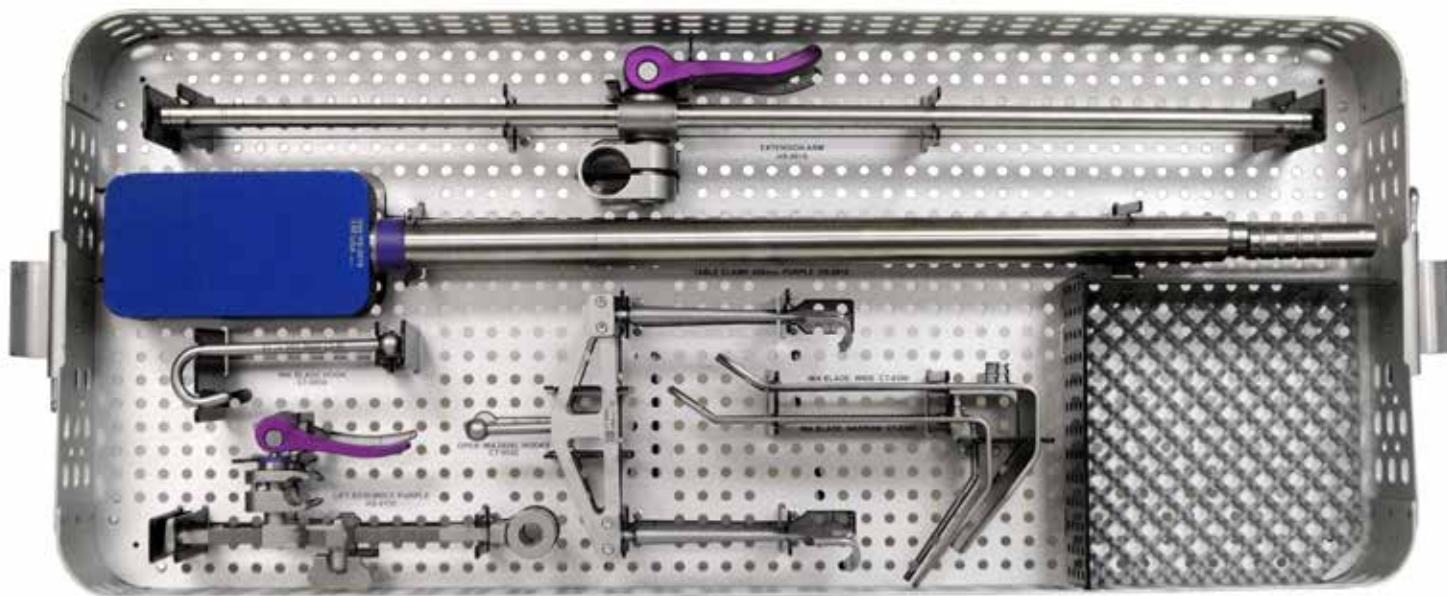

## TeDan Surgical Innovations, LLC

Ph: 713-726-0886 Fax: 713-726-0846 Toll Free: 877-726-0886

12615 West Airport Blvd., Suite #200, Sugar Land, Texas 77478

Email: [info@tedansurgical.com](mailto:info@tedansurgical.com)

[www.tedansurgical.com](http://www.tedansurgical.com)

LIT-0089, REV 3 ©TSI 2017