


CLEAN WAVE KERRISON RONGEUR SET

When performance meets efficiency

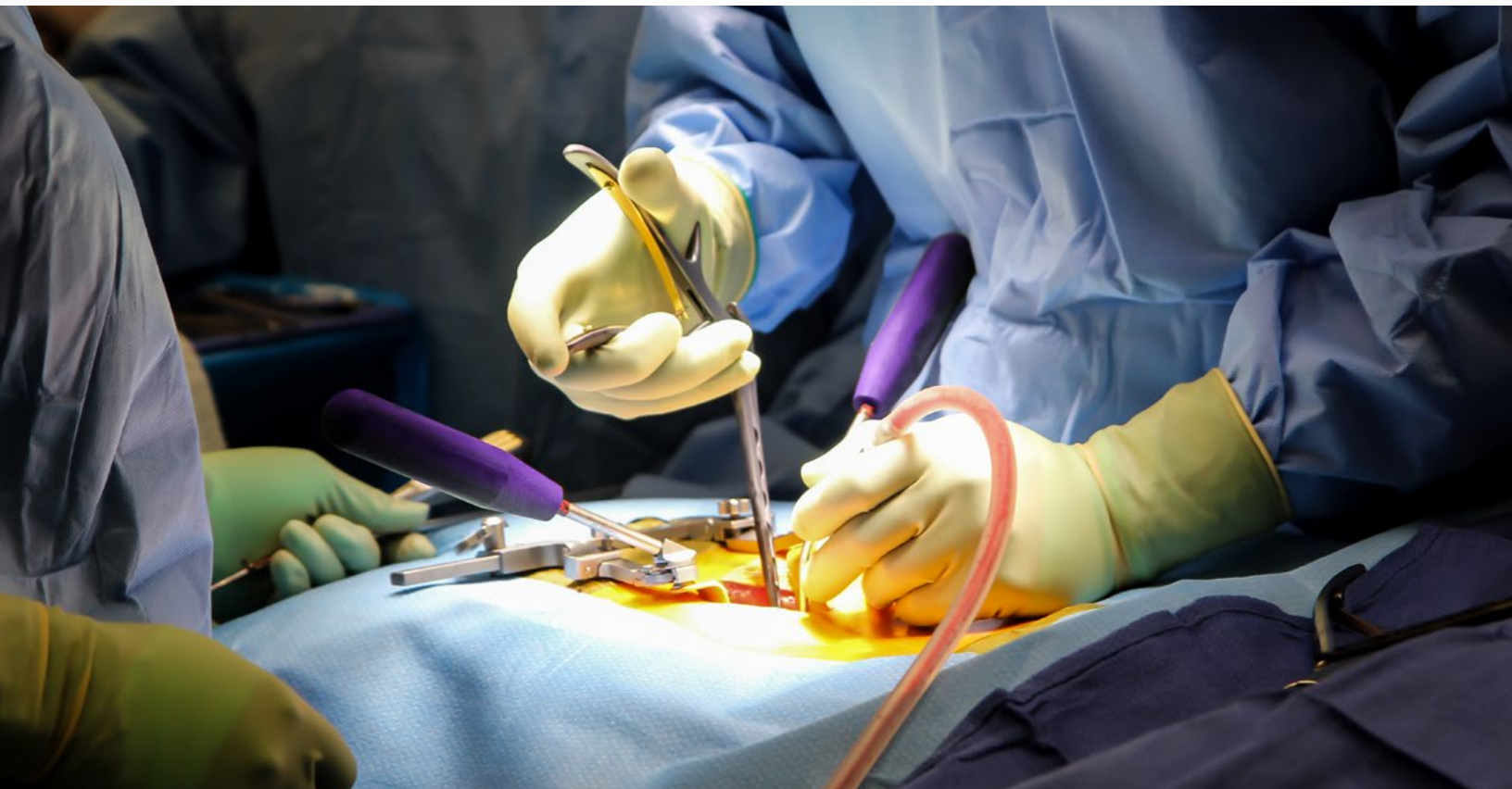

Kerrison Rongeurs, along with other hand held instruments, have long been thought of as some of the most difficult items to reprocess due to their bone collection capabilities and shaft based design. With growing concerns relating to cross contamination, it is more important now than ever to provide easy-to-clean instrumentation.


The Clean Wave Kerrison Rongeur Set from TeDan Surgical Innovations (TSI) delivers the perfect combination of traditional Kerrison bone cutting power with the cleaning efficiency of a take apart instrument—with no disassembly required. With patient safety at the forefront, Clean Wave Kerrison Rongeurs give surgeons the ability to focus on the procedure at hand.

Enabling Effective Cleaning



Distinct wave design increases permeability during the cleaning process, ensuring cross contaminants are removed before each case.

Superior Performance

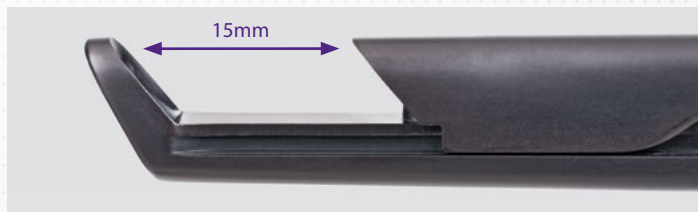

Unique to TSI, a large 15 mm opening* increases working space to provide efficient cutting capability.



Ceramic coating reduces glare and ensures a smooth glide for optimal performance.

CERAMIC COATED CLEAN WAVE KERRISONS

| | |
|----------|--------------------------------------|
| KR-1311B | 8" Kerrison, 1 mm, 15 mm Open, Black |
| KR-1312B | 8" Kerrison, 2 mm, 15 mm Open, Black |
| KR-1313B | 8" Kerrison, 3 mm, 15 mm Open, Black |
| KR-1314B | 8" Kerrison, 4 mm, 15 mm Open, Black |
| KR-1315B | 8" Kerrison, 5 mm, 15 mm Open, Black |
| KR-1316B | 8" Kerrison, 6 mm, 15 mm Open, Black |

STANDARD CLEAN WAVE KERRISONS

| | |
|---------|-------------------------------|
| KR-1311 | 8" Kerrison, 1 mm, 15 mm Open |
| KR-1312 | 8" Kerrison, 2 mm, 15 mm Open |
| KR-1313 | 8" Kerrison, 3 mm, 15 mm Open |
| KR-1314 | 8" Kerrison, 4 mm, 15 mm Open |
| KR-1315 | 8" Kerrison, 5 mm, 15 mm Open |
| KR-1316 | 8" Kerrison, 6 mm, 15 mm Open |

Bayoneted shaft enables maximum visualization for minimally invasive procedures.

CERAMIC COATED BAYONETED CLEAN WAVE KERRISONS 10 MM OPENING

| | |
|----------|--------------------------------------|
| KR-1411B | 9" Bayoneted Kerrison, 1.5 mm, Black |
| KR-1412B | 9" Bayoneted Kerrison, 2 mm, Black |
| KR-1413B | 9" Bayoneted Kerrison, 3 mm, Black |
| KR-1414B | 9" Bayoneted Kerrison, 4 mm, Black |
| KR-1415B | 9" Bayoneted Kerrison, 5 mm, Black |

INSTRUMENT TRAY

| | |
|---------|---------------------------------------|
| KR-0700 | Instrument Tray, Holds 12 Clean Waves |
|---------|---------------------------------------|



Images courtesy of Bonaventure Ngu, MD, F.A.O.S.;
Photo credit: Chris David

*Only applies to 8" Clean Wave Kerrisons



Ph: 713-726-0886 **Fax:** 713-726-0846 **Toll Free:** 877-726-0886
12615 West Airport Blvd., Suite #200, Sugar Land, Texas 77478
Email: info@tedansurgical.com
www.tedansurgical.com
LIT-0073, REV 1 ©TSI 2018